

**NEW
FULL RANGE**

PRODIMED
CRITICAL CARE

SELDIFLEX[®]2
CENTRAL VENOUS CATHETER


SELDIFLEX[®]2 Lumen central venous catheter may be used for:

- Therapeutical agents administration
- Measurement of central venous pressure
- Vascular filling
- Parenteral nutrition

SELDIFLEX[®]2 Multi-lumen central venous catheter may be used for:

- Multiple concomitant medications
- Amines administration
- Antibiotics administration
- Measurement of central venous pressure
- Vascular filling
- Parenteral nutrition
- Hyperosmolar solutions administration (>700 mOsm)

SELDIFLEX[®]2 kits composition

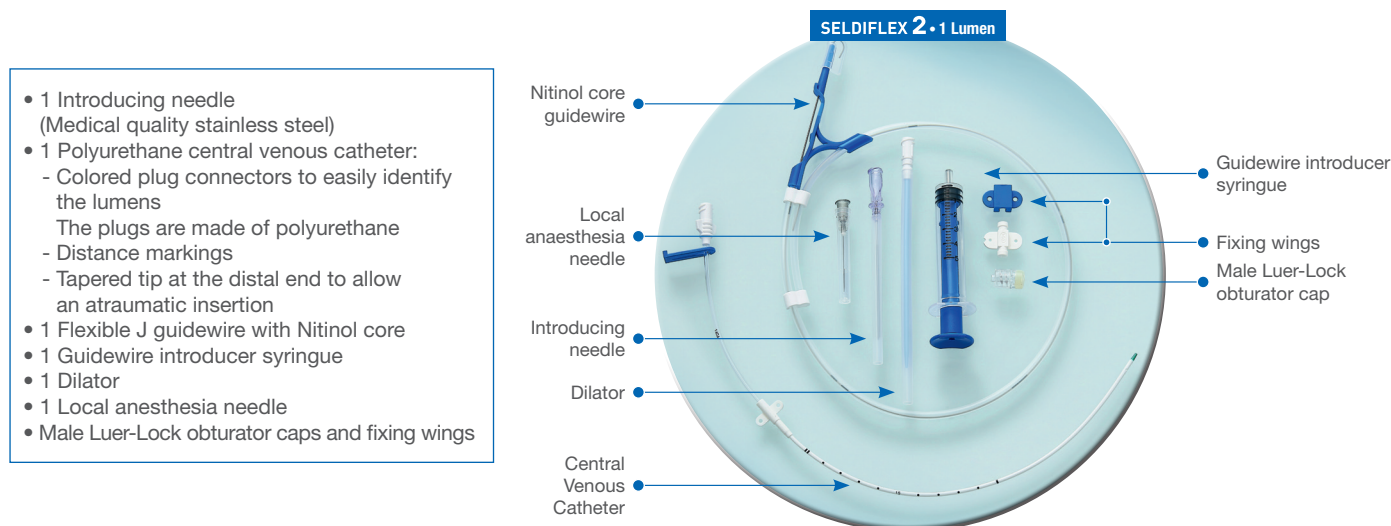

SELDIFLEX[®]2 • 1 Lumen Packaging: 10 / box

Designation	Order Reference	Catheter diameter (French / Gauge)	Catheter effective length (cm)	Dilator size (French x cm)	Needle diameter (Gauge)	Needle length (cm)	Guidewire	Lumen Size (Gauge)	Lumen diameter (Gauge)	Flow Rate (ml/min)	Priming Volume (ml)
Seldiflex 2 1L 5cm 2F	FV-1221	2 Fr (22G)	5 cm	4 Fr x 5 cm	20G	3.3 cm	0.018" x 50 cm	Distal	22G	10-15	0.144
Seldiflex 2 1L 8cm 2F	FV-1222	2 Fr (22G)	8 cm	4 Fr x 5 cm	20G	3.3 cm	0.018" x 50 cm	Distal	22G	8-13	0.151
Seldiflex 2 1L 5cm 3Fr	FV-1321	3 Fr (20G)	5 cm	4 Fr x 5 cm	20G	3.3 cm	0.021" x 50 cm	Distal	20G	32-35	0.16
Seldiflex 2 1L 8cm 3F	FV-1322	3 Fr (20G)	8 cm	4 Fr x 5 cm	20G	3.3 cm	0.021" x 50 cm	Distal	20G	15-30	0.172
Seldiflex 2 1L 10cm 3Fr	FV-1323	3 Fr (20G)	10 cm	4 Fr x 5 cm	20G	3.3 cm	0.021" x 50 cm	Distal	20G	24-27	0.18
Seldiflex 2 1L 16cm 3Fr	FV-1325	3 Fr (20G)	16 cm	4 Fr x 5 cm	20G	3.3 cm	0.021" x 50 cm	Distal	20G	21-23	0.19
Seldiflex 2 1L 20cm 3Fr	FV-1326	3 Fr (20G)	20 cm	4 Fr x 5 cm	20G	3.3 cm	0.021" x 60 cm	Distal	20G	16-21	0.23
Seldiflex 2 1L 16cm 4F	FV-1425	4 Fr (18G)	16 cm	5,5 Fr x 6 cm	19G	3.8 cm	0.025" x 50 cm	Distal	18G	25-35	0.238
Seldiflex2 1L 20cm 4F	FV-1426	4 Fr (18G)	20 cm	5,5 Fr x 6 cm	19G	3.8 cm	0.025" x 60 cm	Distal	18G	20-30	0.252
Seldiflex 2 1L 30cm 4F	FV-1428	4 Fr (18G)	30 cm	5,5 Fr x 6 cm	19G	3.8 cm	0.025" x 70 cm	Distal	18G	15-20	0.321
Seldiflex 2 1L 16cm 5F	FV-1525	5 Fr (16G)	16 cm	6 Fr x 10 cm	18G	6.5 cm	0.032" x 50 cm	Distal	16G	45-65	0.373
Seldiflex 2 1L 20cm 5F	FV-1526	5 Fr (16G)	20 cm	6 Fr x 10 cm	18G	6.5 cm	0.032" x 60 cm	Distal	16G	35-60	0.435
Seldiflex 2 1L 30cm 5F	FV-1528	5 Fr (16G)	30 cm	6 Fr x 10 cm	18G	6.5 cm	0.032" x 70 cm	Distal	16G	30-55	0.498
Seldiflex 2 1L 16cm 6F	FV-1625	6 Fr (14G)	16 cm	8 Fr x 10 cm	18G	6.5 cm	0.032" x 50 cm	Distal	14G	85-120	0.494
Seldiflex 2 1L 20cm 6F	FV-1626	6 Fr (14G)	20 cm	8 Fr x 10 cm	18G	6.5 cm	0.032" x 60 cm	Distal	14G	70-110	0.526

the Seldinger technique

SELDIFLEX[®] 2 • 2 Lumen Packaging: 10 / box

	Designation	Order Reference	Catheter diameter (French / Gauge)	Catheter effective length (cm)	Dilator size (French x cm)	Needle diameter (Gauge)	Needle length (cm)	Guidewire	Lumen Size (Gauge)	Lumen diameter (Gauge)	Flow Rate (ml/ min)	Priming Volume (ml)
	 Seldiflex 2 2L 5cm 4F	FV-2421	4 Fr (18G)	5 cm	5.5 Fr x 6 cm	20G	3.3 cm	0.018" x 50 cm	Distal Proximal	22G 22G	17-22 20-35	0.16 0.17
	Seldiflex 2 2L 8cm 4F	FV-2422	4 Fr (18G)	8 cm	5.5 Fr x 6 cm	20G	3.3 cm	0.018" x 50 cm	Distal Proximal	22G 22G	15-25 9-20	0.16 0.172
	Seldiflex2 2L 13cm 4F	FV-2424	4 Fr (18G)	13 cm	5.5 F x 6 cm	20G	3.3 cm	0.018" x 50 cm	Distal Proximal	22G 22G	12-15 12-15	0.192 0.200
	Seldiflex 2 2L 10cm 5F	FV-2523	5 Fr (16G)	10 cm	6 Fr x 7 cm	20G	3.3 cm	0.021" x 50 cm	Distal Proximal	18G 20G	20-40 15-30	0.181 0.161
	Seldiflex 2 2L 13cm 5F	FV-2524	5 Fr (16G)	13 cm	6 Fr x 7 cm	20G	3.3 cm	0,021" x 50 cm	Distal Proximal	18G 20G	20-35 10-25	0.22 0.201
	Seldiflex2 2L 16cm 5F	FV-2525	5 Fr (16G)	16 cm	6 Fr x 7 cm	20G	3.3 cm	0.021" x 50 cm	Distal Proximal	18G 20G	22-24 22-24	0.21 0.23
	Seldiflex 2 2L 16cm 7F	FV-2725	7 Fr (13G)	16 cm	8 Fr x 10 cm	18G	6.5 cm	0.035 " x 50 cm	Distal Proximal	14G 18G	75-110 25-40	0.472 0.351
	Seldiflex 2 2L 20cm 7F	FV-2726	7 Fr (13G)	20 cm	8 Fr x 10 cm	18G	6.5 cm	0.035 " x 60 cm	Distal Proximal	14G 18G	55-95 15-40	0.505 0.366
	Seldiflex 2 2L 30cm 7F	FV-2728	7 Fr (13G)	30 cm	8 Fr x 10 cm	18G	6.5 cm	0.035 " x 70 cm	Distal Proximal	14G 18G	55-75 10-25	0.65 0.456
	Seldiflex 2 2L 16cm 8.5F	FV-2925	8.5 Fr (12G)	16 cm	9 Fr x 10 cm	18G	6.5 cm	0.035 " x 50 cm	Distal Proximal	14G 14G	75-120 70-110	0.517 0.466
	Seldiflex 2 2L 20cm 8.5F	FV-2926	8.5 Fr (12G)	20 cm	9 Fr x 10 cm	18G	6.5 cm	0.035 " x 60 cm	Distal Proximal	14G 14G	65-110 60-110	0.552 0.526
	Seldiflex 2 2L 30cm 8.5F	FV-2928	8.5 Fr (12G)	30 cm	9 Fr x 10 cm	18G	6.5 cm	0.035 " x 70 cm	Distal Proximal	14G 14G	50-80 50-90	0.692 0.715

SELDIFLEX[®] 2 • 3 Lumen Packaging: 10 / box

	Designation	Order Reference	Catheter diameter (French / Gauge)	Catheter effective length (cm)	Dilator size (French x cm)	Needle diameter (Gauge)	Needle length (cm)	Guidewire	Lumen Size (Gauge)	Lumen diameter (Gauge)	Flow Rate (ml/ min)	Priming Volume (ml)
	 Seldiflex 2 3L 5cm 4.5F	FV-3421	4.5 Fr (17G)	5 cm	6 Fr x 5 cm & 6 Fr x 7 cm	20G	3.3 cm	0.021" x 50 cm	Distal Proximal Médial	23G 23G 20G	29-31 13-15 14-16	0.15 0.18 0.14
	 Seldiflex 2 3L 8cm 4.5F	FV-3422	4.5 Fr (17G)	8 cm	6 Fr x 5 cm & 6 Fr x 7 cm	20G	3.3 cm	0.021" x 50 cm	Distal Proximal Médial	23G 23G 20G	24-31 10-13 9-12	0.18 0.17 0.19
	 Seldiflex 2 3L 10cm 4.5F	FV-3423	4.5 Fr (17G)	10 cm	6 Fr x 5 cm & 6 Fr x 7 cm	20G	3.3 cm	0.021" x 50 cm	Distal Proximal Médial	23G 23G 20G	22-24 9-11 9-11	0.17 0.18 0.18
	Seldiflex 2 3L 8cm 5.5F	FV-3522	5.5 Fr (15G)	8 cm	6,5 Fr x 10cm	20G	3.3 cm	0.021" x 50 cm	Distal Proximal Médial	20G 22G 22G	25-40 15-20 15-20	0.201 0.156 0.18
	 Seldiflex 2 3L 13cm 5.5F	FV-3524	5.5 Fr (15G)	13 cm	6,5 Fr x 10cm	20G	3.3 cm	0.021" x 50 cm	Distal Proximal Médial	22G 22G 20G	26-33 11-14 11-14	0.23 0.23 0.20
	Seldiflex 2 3L 16cm 7F	FV-3725	7F (13G)	16 cm	8 Fr x 10 cm	18G	6.5 cm	0.035" x 50 cm	Distal Proximal Médial	16G 18G 18G	60-80 20-35 20-35	0.392 0.322 0.352
	Seldiflex 2 3L 20cm 7F	FV-3726	7F (13G)	20 cm	8 Fr x 10 cm	18G	6.5 cm	0.035" x 60 cm	Distal Proximal Médial	16G 18G 18G	40-65 13-30 13-30	0.437 0.332 0.345
	Seldiflex 2 3L 30cm 7F	FV-3728	7F (13G)	30 cm	8 Fr x 10 cm	18G	6.5 cm	0.035" x 70 cm	Distal Proximal Médial	16G 18G 18G	35-55 10-20 10-20	0.52 0.37 0.402
	 Seldiflex 2 3L 16cm 8.5F	FV-3925	8.5 Fr (12G)	16 cm	9 Fr x 10 cm	18G	6.5 cm	0.035" x 50 cm	Distal Proximal Médial	16G 16G 14G	54-69 54-69 95-121	0.38 0.42 0.50
	Seldiflex 2 3L 20cm 8.5F	FV-3926	8.5 Fr (12G)	20 cm	9 Fr x 10 cm	18G	6.5 cm	0.035" x 60 cm	Distal Proximal Médial	16G 16G 14G	50-63 50-63 81-104	0.402 0.575 0.465
	 Seldiflex 2 3L 30cm 8.5F	FV-3928	8.5 Fr (12G)	30 cm	9 Fr x 10 cm	18G	6.5 cm	0.035" x 70 cm	Distal Proximal Médial	16G 16G 14G	45-58 45-58 72-92	0.52 0.55 0.72

 New

SELDIFLEX® 2 • 4 Lumen **Packaging: 10 / box**

Designation	Order Reference	Catheter diameter (French / Gauge)	Catheter effective length (cm)	Dilator size (French x cm)	Needle diameter (Gauge)	Needle length (cm)	Guidewire	Lumen Size (Gauge)	Lumen diameter (Gauge)	Flow Rate (ml/min)	Priming Volume (ml)
Seldiflex 2 4L 16 cm 8.5F	FV-4925	8.5 Fr (12G)	16 cm	9 Fr x 10 cm	18G	6.5 cm	0.035" x 50 cm	Distal Proximal Médial 1 Médial 2	16G 18G 14G 18G	40-80 65-105 20-50 20-50	0.350 0.436 0.312 0.329
Seldiflex 2 4L 20cm 8.5F	FV-4926	8.5 Fr (12G)	20 cm	9 Fr x 10 cm	18G	6.5 cm	0.035" x 60 cm	Distal Proximal Médial 1 Médial 2	16G 18G 14G 18G	50-65 65-85 30-35 30-40	0.390 0.420 0.355 0.372
Seldiflex 2 4L 30 cm 8.5F	FV-4928	8.5 Fr (12G)	30 cm	9 Fr x 10 cm	18G	6.5 cm	0.035" x 70 cm	Distal Proximal Médial 1 Médial 2	16G 18G 14G 18G	30-45 45-60 15-25 15-25	0.500 0.517 0.453 0.438




SELDIFLEX® 2 • 2 Lumen



SELDIFLEX® 2 • 3 Lumen



SELDIFLEX® 2 • 4 Lumen

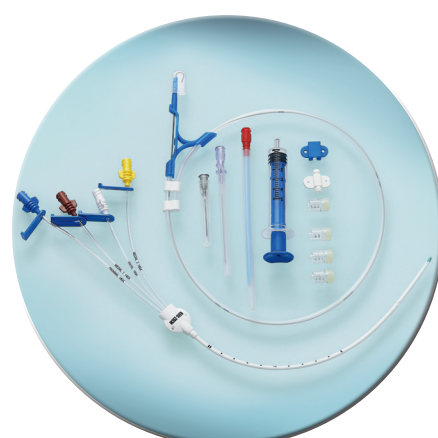

SELDIFLEX® 2 Kits identification


+ EASY TO IDENTIFY

+ EASY

+ EFFICACY

+ SAFETY

Adult version

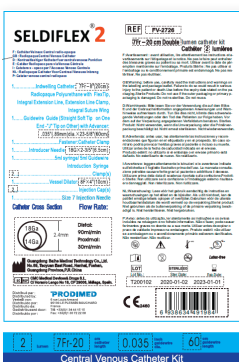


Adult's blue peel-open packaging



Adult's tray (Blue guidewire support)

Pediatric version

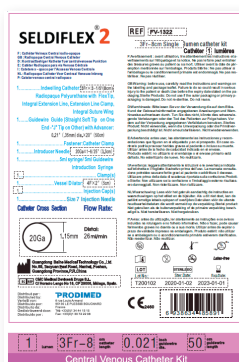

Pediatric's pink peel-open packaging



Pediatric's tray (Translucent guidewire support)

Seldiflex2[®] range

Full range for pediatric and adult patients. **SELDIFLEX[®]2** wide range offers you the right catheter for all kind of procedures. **SELDIFLEX[®]2** allows you to perform the Seldinger technique which is easy and safe to perform your central venous catheterization.

Kit version

SELDIFLEX[®]2 products are always presented in kit. The peel-open packaging colour allows a better differentiation between Pediatric and Adult versions of **SELDIFLEX[®]2**.

All the contents of the kit

- **J guidewire with nitinol core** for a smooth glide and atraumatic passage through the blood vessels.
- **Guidewire support** designed to facilitate gripping even in difficult conditions.
- **Guidewire introducer syringe** for an easy catheterization.
- **Polyuréthane catheter** for a better comfort and a smooth glide. Distance markings allows a better control during the procedure.


Material

SELDIFLEX[®]2 is a polyurethane catheter.

Polyurethane is one of the most polyvalent polymers. It has higher mechanical qualities than silicone (with greater flexibility and higher kink resistance).

High biocompatibility and flexibility for a better efficacy and safety.

Seldinger technique description

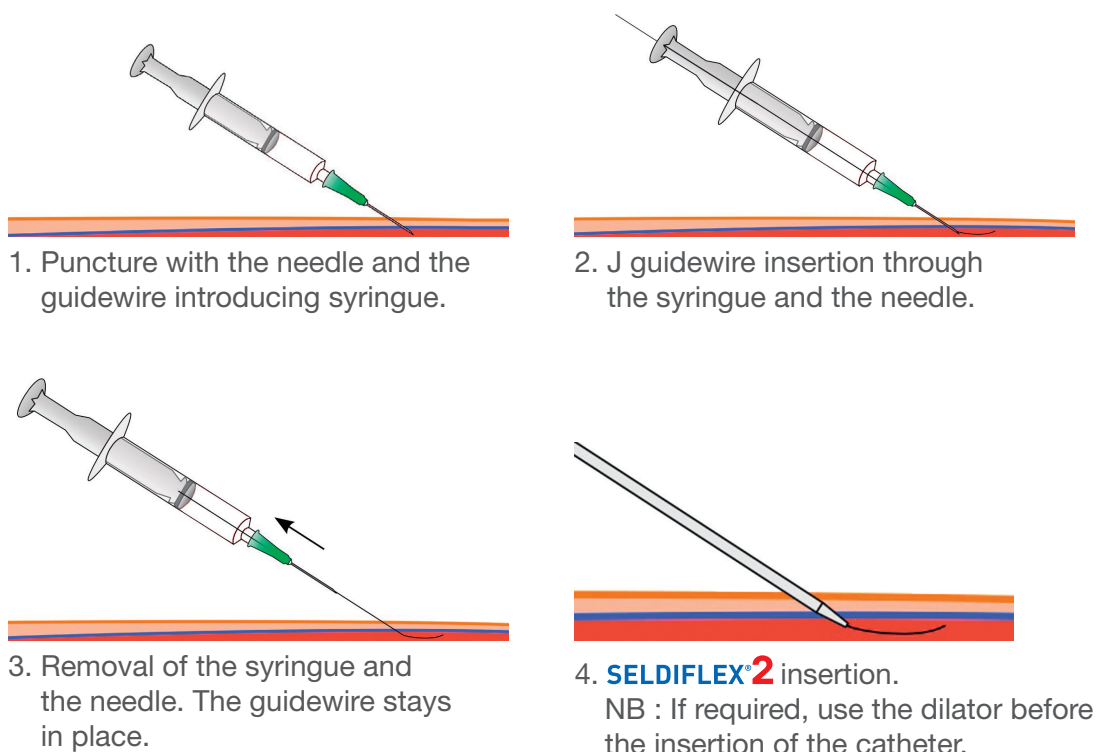

Distributed by:

PRODIMED

6, rue Louis Armand 95130 Le Plessis Bouchard - France

Customer service

Tel. : + 33 (0)1 44 95 17 00 - Fax : +33 (0)1 44 95 14 90

Email : serviceclients@ccdlab.com

www.prodimed.com

Manufactured by:



Guangdong Baihe Medical Technology Co., Ltd.
No.89, Taoyuan East Road, Nanhai, Foshan,
Guangdong Province, P.R.China.



PEELABLE PACKAGING - LATEX FREE - RADIO OPAQUE - SINGLE USE MEDICAL DEVICES

Class III sterile medical devices - CE2460 - Carefully read the instructions. Use before the expiry date indicated on the wrapping. - Do not use if the wrapping or primary packaging is damaged.